

EN

G 1420

Handheld conductivity measuring device



Members of GHM GROUP:

GREISINGER
HONSBERG
Martens
IMTRON
DeltaGHM
VAL.CO

Table of contents

1	About this documentation	4
1.1	Foreword.....	4
1.2	Purpose of the document.....	4
1.3	Legal notices.....	4
1.4	Correctness of content.....	4
1.5	Layout of this document.....	4
1.6	Further information	5
2	Safety	6
2.1	Explanation of safety symbols	6
2.2	Foreseeable misuse	6
2.3	Safety instructions	7
2.4	Intended use	8
2.5	Qualified personnel	8
3	Description	9
3.1	Scope of delivery	9
3.2	Functional description.....	9
4	The product at a glance	10
4.1	The G 1420.....	10
4.2	Display elements	10
4.3	Operating elements	10
5	Bases for measurement	12
5.1	Conductivity principles	12
5.2	Conductivity measurement	12
5.3	Measurement of specific resistance	12
5.4	Electrodes / measuring cell.....	12
5.4.1	Design and selection	12
5.4.2	Calibration / adjustment of the measuring cell.....	13
5.5	Temperature compensation.....	13
5.5.1	NLF temperature compensation according to EN 27888.....	13
5.5.2	Linear temperature compensation LIN	13
6	Operation and maintenance	15
6.1	Operating and maintenance notices	15
6.2	Battery	15
6.2.1	Battery indicator.....	15
6.2.2	Changing battery	15
6.3	Calibration and adjustment service.....	16
6.3.1	Certificates	16
7	Operation	18
7.1	Commissioning	18
7.1.1	Explanation	18
7.2	Configuration	18
7.2.1	Explanation	18
7.2.2	Opening the configuration menu.....	18
7.2.3	Configuring parameters of the configuration menu.....	19
7.2.4	Adjustment of the measuring input	21

7.2.5	Configuring parameters of the adjustment menu.....	21
8	Error and system messages.....	24
9	Disposal.....	25
10	Technical data.....	26
11	Spare parts and accessories.....	27
12	Service.....	28
12.1	Manufacturer.....	28
12.2	Repairs processing.....	28
12.3	Sales offices.....	28
12.4	Sales subsidiaries.....	29

1 About this documentation

1.1 Foreword

Read this document carefully and familiarise yourself with the operation of the product before you use it. Keep this document ready to hand and in the immediate vicinity of the product so that it is available to the personnel/user for reference at all times in case of doubt.

The product was developed according to the state of the art and fulfils the requirements of the applicable European and national Directives. All corresponding documents are available from the manufacturer.

Only technically qualified persons are permitted to carry out installation, commissioning, operation and decommissioning. The qualified personnel must have carefully read and understood the operating manual before beginning any work.

Only technically qualified persons are permitted to carry out commissioning, operation, maintenance and decommissioning. The qualified personnel must have carefully read and understood the operating manual before beginning any work.

1.2 Purpose of the document

- It provides important information for working safely and efficiently with the product.
- In addition to the quick reference guide with all relevant legal and safety content in hard copy, this document is a detailed reference option for the product.

1.3 Legal notices

The liability and warranty of the manufacturer for damages and consequential damages are voided with misuse, disregarding this document, disregarding safety notices, assignment of inadequately qualified technical personnel and arbitrary modifications of the product.

Only carry out the maintenance and service tasks on this product that are described in this documentation. In the process, adhere to the specified steps. For your own safety, only use original spare parts and accessories of the manufacturer. We assume no liability for the use of other products and resulting damage.

This document is entrusted to the recipient for personal use only. Any impermissible transfer, duplication, translation into other languages or excerpts from this operating manual are prohibited.

The manufacturer assumes no liability for print errors.

1.4 Correctness of content

The contents of this document were checked for corrected and are subject to a continuous correction and updating process. This does not rule out potential errors. In the event that errors are discovered or in case of suggestions for improvement, please inform us immediately via the indicated contact information in order to help us make this document even more user-friendly.

1.5 Layout of this document

Description

Each chapter is explained at the beginning in the description.

Prerequisite

All mandatory prerequisites are then listed for each step.

Instruction

Tasks to be carried out by the personnel / user are represented as numbered instructions. Adhere to the sequence of the specified instructions.

Representation

Shows an illustrative instruction or a configuration of the product.

Formula

Some instructions include a formula for a general understanding of a configuration, programming or a setting of the product.

Outcome of an action

Result, consequence or effect of an instruction.

Emphases

In order to simplify legibility and provide a clearer overview, various sections / information are emphasised.

- *Display elements*
- *Mechanical controls*
- **Product functions**
- **Product labels**
- Cross-reference [▶ p. 4]
- *Foot notes*

1.6 Further information

Software version of the product:

- V1.3 or later

2 Safety

2.1 Explanation of safety symbols



DANGER

This symbol warns of imminent danger which can result in death, severe bodily injury, or severe property damage in case of non-observance.



CAUTION

This symbol warns of potential dangers or harmful situations which can cause damage to the device or to the environment in case of non-observance.



CAUTION

This symbol warns of electrical voltage.



NOTE

This symbol indicates processes which can have a direct influence on operation or can trigger an unforeseen reaction in case of non-observance.

2.2 Foreseeable misuse

The fault-free function and operational safety of the product can only be guaranteed if generally applicable safety precautions and the device-specific safety instructions for this document are observed.

If these notices are disregarded, personal injury or death, as well as property damage can occur.



DANGER

Incorrect area of application!

In order to prevent erratic behaviour of the product, personal injury or property damage, the product must be used exclusively as described in the chapter Description [► p. 9] in the operating manual.

- Do not use in safety / Emergency Stop devices!
- The product is not suitable for use in explosion-prone areas!
- The product must not be used for diagnostic or other medical purposes on patients!
- The product is not intended to come into direct contact with food. For measurement in foods, samples must be taken and discarded after the measurement!
- Not suitable for use with requirements on functional safety, e.g. SIL!



CAUTION

Break-out

In order to prevent the mount or product itself from breaking out, stress must be avoided.

- Do not hang anything on it!
- Do not place anything on it!

2.3 Safety instructions

This product has been designed and tested according to the safety requirements for electronic measuring devices. The product must be used according to the technical data. Technical data [▶ p. 26].



CAUTION

Impaired function

Errors during installation, assembly or configuration of GHM products can result in impaired function or damage in the downstream process. Severe danger, however, does not emanate directly from our products.

- Ensure proper and correct installation, assembly and programming and configuration of the product!
- Provide independent safety devices!
- Settings must be made by qualified personnel only!



CAUTION

Current strength

The maximum current strength, which is specified in the chapter Technical data [▶ p. 26], must not be exceeded.

- Use the correct mains adapter!
- Ensure that the wiring and connections are correct!



CAUTION

Erratic behaviour!

On suspicion that the product can no longer be operated without danger, it must be decommissioned and prevented from recommissioning with appropriate labelling. The safety of the user can be impaired by the device if, for example, if it shows visible damage, it no longer works as specified or if it was stored for an extended period of time under unsuitable conditions.

- Visual inspection!
- In case of doubt, send the product to the manufacturer for repair or maintenance!



CAUTION

Connection error!

A defect can occur if the product is installed without consideration of the thread or connection type.

- The connection flange/thread must be correct!
- Ensure that the mating plug has the correct design!
- Use of capacitors to prevent expanded pipelines due to heat!



NOTE

This product does not belong in children's hands!



NOTE

Installation in zones

Device versions with an ATEX approval are approved for use as an operating device for connection of intrinsically safe sensors installed in Zones 0 or 1, 20 or 21. In the process, the technical safety data must be observed.

2.4 Intended use

The product is designed for measuring the conductivity in liquids. The measuring cell is connected permanently.

2.5 Qualified personnel

For commissioning, operation and maintenance, the relevant personnel must have adequate knowledge of the measuring process and the significance of the measurements. This document makes a valuable contribution to this. The instructions in this document must be understood, observed and followed.

In order to avoid any risks arising from interpretation of the measurements in the concrete application, the user must have additional expertise. The user is solely liable for damages/danger resulting from misinterpretation due to inadequate expertise.

3 Description

3.1 Scope of delivery

Please check to ensure the completeness of the product after opening the package. You should find the following components:

- Quick reference guide
- Handheld measuring device, ready for operation, including batteries
- Permanently connected conductivity measuring cell
- Test report

3.2 Functional description

The product offers precision, speed and reliability in a compact, ergonomic housing. Additional impressive features include the dust-proof and waterproof design in accordance with IP 65/67 and the 3-line illuminated display, which offers overhead display at the push of a button. The product can be switched on, switched off and configured and the measurements and parameters can be adjusted and held with the operating elements. The product is equipped with a universal 2-pole stainless steel conductivity measuring cell optimised for purest for measurement in the range of 0.000 $\mu\text{S}/\text{cm}$ to 100.0 $\mu\text{S}/\text{cm}$. In addition to conductivity, specific resistance can be determined.

4 The product at a glance

4.1 The G 1420



4.2 Display elements

Display



Battery indicator

Evaluation of the battery status



Unit display

Display of units, if applicable, with unstable symbol or type of mode, min/max/hold



Main display

Measurement of the current conductivity value or value for min/max/hold



Auxiliary display

Corresponding temperature value for the value shown in the main display. If applicable, alternating with the temperature compensation.



NOTE

The unit display shows a rotating circle segment in the first position as long as the measurement is unstable, if the position is unoccupied by the unit display.

4.3 Operating elements



On / Off button

Press briefly

Switch on the product

Activate / deactivate lighting

Long press

Switch off the product

Reject changes in a menu

**Up / Down button**

Press briefly

Display of the min/max value

Long press

Change value of the selected parameter

Both simultaneously

Reset the min/max value of the current measurement

Rotate display, overhead display

**Function key**

Press briefly

Freeze measurement

Return to measurement display

Call up next parameter

Long press, 2s

Open menu, frozen measurement is displayed

Close menu, changes are saved

5 Bases for measurement

5.1 Conductivity principles

Conductivity γ

Conductivity describes the capability of a material to conduct electrical current. It is also the inverse of specific resistance. The conductance is the inverse of the measured resistance R.

Formula

$$\gamma = l / (R \cdot A)$$

l = length of the material

A = cross section

R = measured resistance

Unit [γ] = Siemens / metre = S / m

Normally, the values for liquids are specified in $\mu\text{S} / \text{cm}$ or in mS / cm . With measurements in purest water, the units $\text{k}\Omega \cdot \text{cm}$ or $\text{M}\Omega \cdot \text{cm}$ are commonly used.

5.2 Conductivity measurement

The conductivity measurement is a comparatively uncomplicated measurement. The standard electrodes are stable for correct use for a long time and can be calibrated over the gradient correction.

Range	1	2	3
	0.000 to 2.000 $\mu\text{S} / \text{cm}$	0.00 to 20.00 $\mu\text{S} / \text{cm}$	0.0 to 100.0 $\mu\text{S} / \text{cm}$

With the automatic range selection, the range with the best resolution is selected automatically.

5.3 Measurement of specific resistance

The specific resistance is the inverse of the conductivity and is specified in the device in $\text{M}\Omega \cdot \text{cm}$ or $\text{k}\Omega \cdot \text{cm}$.

Range	1	2	3
	10.0 to 200.0 $\text{k}\Omega \cdot \text{cm}$	0.010 to 2.000 $\text{M}\Omega \cdot \text{cm}$	0.01 to 20.00 $\text{M}\Omega \cdot \text{cm}$


With the automatic range selection, the range with the best resolution is selected automatically.

5.4 Electrodes / measuring cell

5.4.1 Design and selection

There are basically two different types of measuring cells: 2-pole and 4-pole measuring cells. Control and/or evaluation take place in a similar manner; the 4-pole measuring cells can compensate well for polarisation effects and contamination to a certain degree with the more elaborate measuring processes.

The product is equipped with a permanently connected 2-pole measuring cell.

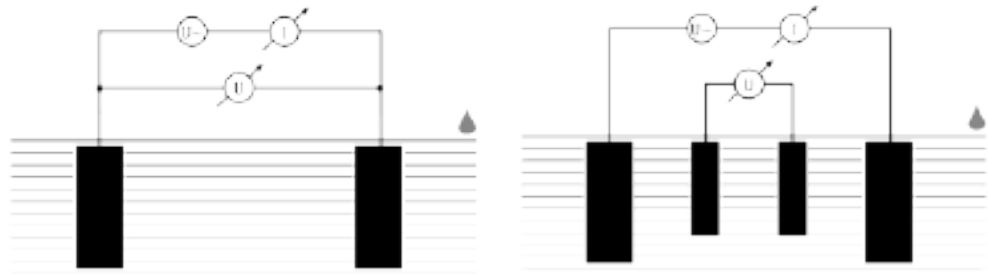

2-pole measuring cell

4-pole measuring cell

2-pole stainless steel measuring cell for low conductivities $<100\mu\text{S}$, e.g. pure and purest water, boiler water, osmosis and filter technology

5.4.2 Calibration / adjustment of the measuring cell

In harsh applications and due to ageing processes, the cell constant of measuring cells changes. Depending on the application and precision requirement, the overall precision of the display device and measuring cell measuring chain should be checked regularly. Special testing and calibration solutions, such as GKL 101 are available for this purpose. In normal application conditions, semi-annual testing is recommended; Adjustment of the measuring input [► p. 21]. A system test by the manufacturer is recommended in case of doubt; Calibration and adjustment service [► p. 16].

For this purpose, also refer to

- ▣ Adjustment of the measuring input [► 21]
- ▣ Calibration and adjustment service [► 16]

5.5 Temperature compensation

The conductivity of aqueous solutions is temperature-dependent. The temperature dependency varies greatly according to the type of solution. With temperature compensation, the solution is calculated back to a uniform temperature in order to compare it independently of the temperature. The normal operating temperature for this is 25 °C. However, it can also be adjusted to 20 °C.

5.5.1 NLF temperature compensation according to EN 27888

For most applications, such as fish husbandry applications and measurement of surface water and drinking water, non-linear temperature compensation for natural water nLF is sufficiently accurate in accordance with EN 27888.

The normal operating temperature is 25 °C.

The recommended application range of nLF compensation is between 60 $\mu\text{S}/\text{cm}$ and 1000 $\mu\text{S}/\text{cm}$.

5.5.2 Linear temperature compensation LIN

If the function of the temperature compensation is not known exactly, a linear temperature compensation is normally adjusted in the product.

Open the **Configuration** menu and select the parameter $tcor$. Configuration parameters LIn and $ELIn$ correspond to TK_{in} .

Put simply, this means that the temperature dependency is approximately the same over the considered concentration range of the solution.

Temperature coefficients of 2.0 %/K are most common.

Formula

$$LF_{T_{ref}} = LF_{TX} / ((1 + TK_{lin} / 100\%) * (TX - T_{ref}))$$

A temperature coefficient can be determined, for example, with measurement of a solution at 2 temperatures, T1 and T2, with temperature compensation switched off.

Formula

$$TK_{lin} = ((LF_{T1} - LF_{T2}) * 100\%) / ((T1 - T2) * LF_{T1})$$

TK_{lin} = value is entered in the **Configuration** menu in parameter $tLin$.

LF_{T1} = conductivity at Temperature 1

LF_{T2} = conductivity at Temperature 2

6 Operation and maintenance

6.1 Operating and maintenance notices



NOTE

The product and conductivity measuring cell must be handled with care and used in accordance with the technical data. Do not throw or strike.



NOTE

If the product is stored at a temperature above 50 °C, or is not used for an extended period of time, the batteries must be removed. Leaks from the batteries are avoided as a result.

The device is calibrated at the factory with the permanently connected conductivity measuring cell. The highest system precision can be achieved in this manner. If desired, a gradient correction can be carried out for the product in order to further optimise the accuracy in a narrow range. This is only necessary for normal use. See Adjustment of the measuring input [► p. 21].

6.2 Battery

6.2.1 Battery indicator

If the empty frame in the battery display blinks, the batteries are depleted and must be replaced. However, the device will still operate for a certain length of time.

If the **bAt** display text appears in the main display, the battery voltage is no longer adequate for operation of the product. The battery is fully depleted.

6.2.2 Changing battery



DANGER

Danger of explosion!

Using damaged or unsuitable batteries can generate heat, which can cause the batteries to crack and possibly explode!

- Only use high-quality and suitable alkaline batteries!



CAUTION

Damage!

If the batteries have different charge levels, leaks and thus damage to the product can occur.

- Use new, high-quality batteries!
- Do not use different types of batteries!
- Remove depleted batteries and dispose of them at a suitable collection point!



NOTE

Unnecessary screwing places the water-tightness of the product, among other things, at risk and should be avoided.



NOTE

Read the following handling instructions before replacing batteries and follow them step by step. If disregarded, the product could be damaged or the protection from moisture could be diminished.

Description

Prerequisites

Instruction

Proceed as follows to replace the batteries.

- The product is switched off.

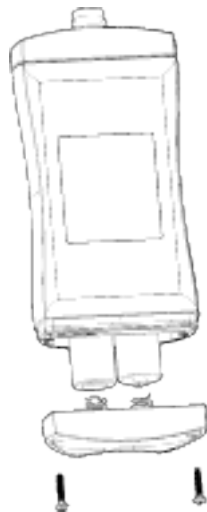


Fig. 4: Changing battery

1. Unscrews the Phillips screws and remove the cover.
2. Carefully replace the two Mignon AA batteries. Ensure that the polarity is correct! It must be possible to insert the batteries in the correct position without using force.
3. The O-ring must be undamaged, clean and positioned at the intended depth. In order to facilitate assembly and avoid damage, a suitable grease can be applied.
4. Fit the cover on evenly. The O-ring must remain at the intended depth!
5. Tighten the Phillips screws.

Outcome of an action

The product is now ready for use again.

6.3 Calibration and adjustment service

6.3.1 Certificates

The certificates are categorised as ISO calibration certificates and DAkkS calibration certificates. The purpose of the calibration is to verify the precision of the measuring device by comparing it with a traceable reference.



NOTE

The ISO standard 9001 is applied for the calibration certificates. These certificates are an affordable alternative to the DAkkS calibration certificates and provide information of the traceable reference, a list of individual values and documentation.



NOTE

The DAkkS calibration is based on DIN EN ISO/17025, the accreditation basis recognised worldwide. These certificates offer high-quality calibration and consistently high quality. DAkkS calibration certificates can only be issued by accredited calibration

laboratories which have demonstrated their expertise in accordance with DIN EN ISO/IEC 17025. The DAkkS calibration includes any necessary adjustment with the purpose of minimising a deviation of the measuring device.

DAkkS calibration certificates are accompanied with a list of individual measurements before and after the adjustment, documentation and, if applicable, graphic representation, calculation of the expanded measuring uncertainty and traceability to the national standard.



NOTE

The product is delivered with a test report. This confirms that the measuring device has been adjusted and tested.



NOTE

Only the manufacturer can check the basic settings and make corrections if necessary.

7 Operation

7.1 Commissioning

7.1.1 Explanation

Description The product is switched on with the *On/Off button*. It may be necessary to configure the product after switching on. See Configuration [► p. 18].

Prerequisite – Sufficiently full batteries are inserted in the product.

Instruction – Press the *On/Off button*.

Outcome of an action Information about the configuration of the product appears in the display.

P_{OFF}	Automatic shut-off	Automatic shut-off activated. The product is switched off if no buttons have been pressed after the adjusted time
ϵ_{oF}	Zero point correction	If a zero point correction of the temperature sensor was made
ϵ_{5L}	Gradient correction	If a gradient correction of the temperature sensor was made
ϵ_{CL}	Gradient correction	If a gradient correction of the conductivity measuring cell was made

– The product is now ready for measurement.



NOTE

The device is calibrated at the factory and ready for measurement. The offset and gradient correction of the temperature measurement and gradient correction of the conductivity measurement can be carried out when exact references are available. This is only necessary in exceptional cases.

7.2 Configuration

7.2.1 Explanation

The following steps describe how to adapt the product for your purposes.



NOTE

There are various configuration parameters available depending on the product version and configuration. They can differ depending on the product version and configuration.

7.2.2 Opening the configuration menu

Description In order to configure the product, you must first open the **Configuration** menu. The menu is opened as shown in the illustration.






Prerequisite 1. Press the *Function key* for 2 seconds to open the **Configuration** menu.

Instruction 2. ϵ_{oF} appears in the display. Release the function key.

3. By briefly pressing the *Function key*, you can scroll through the parameters. Select the parameter you would like to configure.

4. When you have selected the desired parameter, change the parameter to the desired value with the *Up button* and the *Down button*.
5. The changes are saved after running through the entire **Configuration** menu. *Stop* appears in the display. The **Configuration** menu can be exited from any arbitrary parameter by pressing and holding the *Function key* for 2 seconds. The changes made up that point are saved.

Representation

Call up menu	Next parameter	Change value	Save changes	Discard changes
				
2s		Press: Single step Hold: Rapid change	2s	2s Product is switched off

Outcome of an action

The **Configuration** menu is closed after the last parameter.



NOTE

If the product is switched off without saving the configuration, the last save value is reproduced on the next start-up of the product.

7.2.3 Configuring parameters of the configuration menu

Description

The following representation shows the available parameters and various configuration options.




Prerequisite

- The **Configuration** menu is open. See Opening the configuration menu [► p. 18].

Instruction

1. Select the desired parameter you would like to configure.
2. Adjust the desired configuration in the selected parameter with the *Up button* and *Down button*.
3. The available configuration options are listed for each parameter in the following representation.

Representation

Parameter	Values	Meaning
	 	
Input		
<i>inP</i>	<i>Cond</i>	Measured variable - conductivity
	<i>rES</i>	Measured variable - specific resistance
Temperature compensation		
<i>tcor</i>	<i>oFF</i>	Do not compensate conductivity measurement
	<i>nLF</i>	Non-linear function for natural water in accordance with EN 27888 (ISO 7888) Ground water, surface water or drinking water
	<i>nRcl</i>	Compensation of weak NaCl solutions only in pure and ultrapure water
	<i>Li n</i>	Linear temperature compensation

 Compensation coefficient

$\epsilon_{L,n}$	0.300 .. 3.000	Temperature compensation coefficient in %/K
------------------	----------------	---

Reference temperature for temperature compensation

t_{REF}	25 °C	Reference temperature 25 °C or 77 °F
	20 °C	Reference temperature 20 °C or 68 °F

Shut-off time

P_{OFF}	oFF	No automatic shut-off
	15 30 60 120 240	Automatic shut-off after a selected time in minutes, during which no buttons have been pressed

Backlighting

L, tE	oFF	Backlighting deactivated
	15 30 60 120 240	Automatic shut-off of the backlighting after a selected time in seconds, during which no buttons have been pressed
	oN	No automatic shut-off of the backlighting

Temperature unit

$Unit$	°C	Temperature display in °C
	°F	Temperature display in °F

Factory settings

$init$	no	Use current configuration
	YES	Reset product to factory settings. <i>init done</i> appears in the display

Outcome of an action

The changed value is saved and the **Configuration** menu is closed. *Star* appears in the display. If necessary, the product is restarted automatically in order to adopt the changed values.


NOTE

The configuration is closed if no button is pressed for 2 minutes. Any changes made up to that point are not saved. *end* appears in the display.

7.2.4 Adjustment of the measuring input

Description The temperature input can be adjusted with the zero point correction and the gradient correction. If an adjustment is made, you change the pre-adjusted factory settings. This is signalled with the *L.O.F*, *L.S.L* or *S.C.L* display text when the product is switched on. The standard settings of the zero point value and the gradient value of the temperature input is *0.00*. The standard setting of the gradient value of the conductivity value is *1.000*. It signals that no correction is made.

In order to adjust the product, you must first open the **Adjustment** menu. The menu is opened as shown in the illustration.

- Prerequisites**
- Sufficiently full batteries are inserted in the product.
 - The product is switched off.
 - Ice water, regulated precision water baths or a water bath with a reference measurement are available as a reference.

- Instruction**
1. Press and hold the *Down button*.
 2. Press the *On/Off button* to switch on the product and open the **Configuration** menu. Release the *Down button*. The display shows the first parameter.
 3. By briefly pressing the *Function key*, you can scroll through the parameters. Select the parameter you would like to configure.
 4. When you have selected the desired parameter, change the parameter to the desired value with the *Up button* and the *Down button*.
 5. In order to save the new parameter value, press and hold the *Function key* for longer than 1 second.

Representation

Call up menu



Hold



Release

Outcome of an action

The **Configuration** menu is closed after the last parameter.



NOTE

If the product is switched off without saving the configuration, the last save value is reproduced on the next start-up of the product.

7.2.5 Configuring parameters of the adjustment menu

Description The following representation shows the available parameters and various configuration options.

Prerequisites The **Adjustment** menu is open. See Adjustment of the measuring input [▶ p. 21].

- Instruction**
1. Select the desired parameter you would like to configure.
 2. Adjust the desired configuration in the selected parameter with the *Up button* and *Down button*.
 3. The available configuration options are listed for each parameter in the following representation.

Representation

Parameter	Values	Meaning

Zero point correction

 ϵ_{ZF}

0.00

No zero point correction

-5.00 .. 5.00

Zero point correction in °C. and/or at °F -9.00 .. 9.00

Gradient correction of the temperature

 ϵ_{5L}

0.00

No gradient correction

-5.00 .. 5.00

Gradient correction in %

Gradient correction for the conductivity value

 ϵ_{5L}

1.000

No gradient correction

0.800 .. 1.200

Multiplier for the gradient correction

Formula

Zero point correction:

Displayed value = measured value – ϵ_{ZF}

Gradient correction °C:

Display = (measured value – ϵ_{ZF}) * (1 + ϵ_{5L} / 100)

Gradient correction °F:

Display = (measured value – 32 °F – ϵ_{ZF}) * (1 + ϵ_{5L} / 100) + 32 °FGradient correction γ :Display = measured value / ϵ_{5L}

Example calculation

Temperature

- Zero point correction ϵ_{ZF} to 0.00
- Gradient correction ϵ_{5L} to 0.00
- Display unit *Unit* to °C
- Display in ice water -0.2 °C
- Display in ice water setpoint ϵ_{ZF} = 0.0 °C
- Display in water bath 36.6 °C
- Display in water bath setpoint ϵ_{5L} = 37.0 °C
- ϵ_{ZF} = display zero point correction – setpoint zero point
- ϵ_{ZF} = -0.2 °C – 0.0 °C = -0.2 °C
- ϵ_{5L} = (setpoint gradient correction / (display gradient correction – ϵ_{ZF}) – 1) * 100
- ϵ_{5L} = (37.0 °C / (36.6 °C – (-0.2)) - 1) * 100 = 0.54

Conductivity

- Gradient correction ϵ_{5L} to 1.000
- Temperature compensation ϵ_{corr} to OFF
- Control solution GKL 101 as reference
- Display in GKL 101 at 25 °C, setpoint = 84 $\mu\text{S}/\text{cm}$
- Display = 82.5 $\mu\text{S}/\text{cm}$
- ϵ_{5L} = setpoint / display value
- ϵ_{5L} = 84 $\mu\text{S}/\text{cm}$ / 82.5 $\mu\text{S}/\text{cm}$ = 1.018

**NOTE**

The most precise results can be achieved when the control solution temperature is adjusted to 25 °C. With different temperatures, the table value of the corresponding temperature must be used as a setpoint.

Outcome of an action

The changed value is saved and the *Configuration* menu is closed.

**NOTE**

If the product is switched off without saving the configuration, the last save value is reproduced on the next start-up of the product.

8 Error and system messages

Display	Meaning	Possible causes	Remedy
----	Range switching or measured value unstable Measurement far outside of the measuring range	Measuring cell defect Contamination or air bubbles	Wait for the transient effect of the controller Measurement leaves the permissible range Send in for repair
5Err5 Error	No probe or sensor connected Sensor or probe defect Measuring range exceeded or undercut	Sensor or probe missing Defective sensor or probe Incorrect sensor type selected	Connect sensor or probe Connect different sensor or probe Readjust measuring range
No display, unclear characters or no response when buttons are pressed	Battery depleted System error Product is defective	Battery depleted Error in the product Product is defective	Replace battery Send in for repair
bAt	Battery depleted	Battery depleted	Replace battery
Err.1	Measuring range exceeded	Measurement too high Measuring cell defect	The measurement is above the permissible range Check the measuring cell Send in for repair
Err.2	Measuring range is undercut	Measurement too low Measuring cell defect	The measurement is below the permissible range Check the measuring cell Send in for repair
555 Err	System error	Error in the product	Switch product on/off Replace batteries Send in for repair

9 Disposal

Separation by material and recycling of device components and packaging must take place at the time of disposal. The valid regional statutory regulations and directives applicable at the time must be observed.



NOTE

The device must not be disposed of with household waste. Return it to us, freight pre-paid. We will then arrange for the proper and environmentally-friendly disposal.

Private end users in Germany have the possibility of dropping off the product at the municipal collection centre.

Please dispose of empty batteries at the collection points intended for this purpose.



NOTE

Fill in the return form available from the information base online at www.ghm-group.de and sent it in with the product.

10 Technical data

Measuring range	Conductivity	0.000 .. 2.000 $\mu\text{S/cm}$ 0.00 .. 20.00 $\mu\text{S/cm}$ 0.0 .. 100.0 S/cm
	Specific resistance	10.0 .. 200.0 $\text{k}\Omega/\text{cm}$ 0.010 .. 2.000 $\text{M}\Omega/\text{cm}$ 0.01 .. 20.00 $\text{M}\Omega/\text{cm}$
	Salinity	-
	TDS	-
	Temperature	-5.0 .. +105.0 °C (23.0 .. +221.0 °F) – the conductivity measuring cells can be exposed temporarily to temperatures of up to 100 °C and permanently to temperatures of up to 80 °C.
Accuracy	Conductivity	Typ. $\pm 1\%$ of measured value $\pm 0.5\%$ FS
	Temperature	$\pm 0.3\text{ °C}$
Measuring cycle		approx. 10 measurements per second Updating of the display approx. 2 times per second
Display		3-line segment LCD, additional symbols, illuminated (adjustable white, permanent illumination)
Additional functions		Min/Max/Hold
Compensation		Offset and gradient correction - temperature, Gradient correction - conductivity
Housing		Break-proof ABS housing
	Protection rating	IP65 / IP67
	Dimensions L*W*H [mm] and weight	108 * 54 * 28 mm without measuring cell or kink protection 210 g, incl. battery and measuring cell
Operating conditions		-20 to 50 °C; 0 to 95 % r.h. (temporarily 100 % r.h.)
Storage temperature		-20 to 70 °C
Power supply		2*AA battery (included in the scope of delivery)
	Current requirement/ battery life	approx. 2.2 mA, approx. 3.5 mA with lighting Service life > 1000 hours with alkaline batteries (without back-lighting)
	Battery indicator	4-stage battery status indicator, Replacement indicator for depleted batteries: "BAT"
Auto-power-OFF function		The device switches off automatically if this is activated
Directives and standards		The devices conform to the following Directives of the Council for the harmonisation of legal regulations of the Member States: 2014/30/EU EMC Directive 2011/65/EU RoHS Applied harmonised standards: EN 61326-1:2013 Emission limits: Class B Immunity according to Table 2 Additional errors: < 1 % FS EN 50581:2012 The device is intended for mobile use and/or stationary operation in the scope of the specified operating conditions without further limitations.

11 Spare parts and accessories

A selection of spare parts and accessories for this product is listed below.

Article

Number	Name	Description
610049	Mignon battery AA	Mignon AA spare battery
603499	GWZ-01	Flow-through vessel (for measuring cells with Ø 12 mm and Ø 6 mm hose connection)
700040	HD-22-3	Freely positionable, flexible laboratory electrode holding arm for probes with Ø 12 mm
610049	GB AA	AA battery
611373	ST-G1000	Device protection bag with 1 round cut-out

A complete list of all accessories and spare parts is available in our product catalogue or on our home page. We can also provide further information by phone.

Contact

Internet: www.greisinger.de

Tel: +49 94029383-52

12 Service

12.1 Manufacturer

If you have any questions, please do not hesitate to contact us:

Contact

GHM Messtechnik GmbH
GHM GROUP - Greisinger
Hans-Sachs-Str. 26
93128 Regenstauf | GERMANY
Email: info@greisinger.de | www.greisinger.de
WEEE reg. no. DE 93889386



12.2 Repairs processing

Defective products are repaired professionally and quickly in our service centre.

Open hours and contact

Monday to Thursday from 8:00 to 16:00
Friday from 8:00 to 13:00
GHM Messtechnik GmbH
GHM GROUP - Greisinger
Hans-Sachs-Str.26
Service Centre
93128 Regenstauf | GERMANY
Tel: +49 94029383-39
Fax: +49 94029383-33
service@greisinger.de



NOTE

Fill in the return form available from the information base online at www.ghm-group.de and sent it in with the product.

12.3 Sales offices

North Sales Office

Post code: 00000 – 25999 | 27000 – 34999
37000 – 39999 | 98000 – 99999
Email: vertrieb-nord@ghm-messtechnik.de
Tel: +49 4067073-0
Fax: +49 4067073-288

West Sales Office

Post code: 26000 – 26999 | 35000 – 36999
40000 – 69999
Email: vertrieb-west@ghm-messtechnik.de
Tel: +49 2191 9672-0
Fax: +49 2191 9672-40

South Sales Office

Post code: 70000 – 97999
Email: vertrieb-sued@ghm-messtechnik.de
Tel: +49 9402 9383-52
Fax: +49 9402 9383-33

12.4 Sales subsidiaries

Austria
GHM Messtechnik GmbH
Office Austria
Breitenseer Str. 76/1/36
1140 Vienna | AUSTRIA
Phone +43 660 7335603
a.froestl@ghm-messtechnik.de

Brazil & Latin America
GHM Messtechnik do Brasil Ltda
Av. José de Souza Campos, 1073, cj 06
Campinas, SP
13025 320 | BRAZIL
Phone +55 19 3304 3408
Info@grupoghm.com.br

Czech Republic / Slovakia
GHM Greisinger s.r.o.
Ovci hajek 2 / 2153
158 00 Prague 5
Nove Butovice | CZECH REPUBLIC
Phone +420 251 613828
Fax +420 251 612607
info@greisinger.cz | www.greisinger.cz

Denmark
GHM Maaleteknik ApS
Maarslet Byvej 2
8320 Maarslet | DENMARK
Phone +45 646492- 00
Fax +45 646492- 01
info@ghm.dk | www.ghm.dk

France
GHM GROUP France SAS
Parc des Pivolles
9 Rue de Catalogne
69150 Décines-Charpieu (Lyon) | FRANCE
Phone +33 4 72 37 45 30
a.jouanilou@ghm-group.fr

India
GHM Messtechnik India Pvt Ltd.
209 | Udyog Bhavan | Sonowala Road
Gregaon (E) | Mumbai - 400 063
INDIA
Phone +91 22 40236235
info@ghmgroup.in | www.ghmgroup.in

Italy for Greisinger & Delta OHM
GHM GROUP – Delta OHM
Via Marconi 5
35030 Caselle di Selvazzano
Padova (PD) | ITALY
Phone +39 049 8977150
a.casati@ghm-messtechnik.de

Italy for Honsberg, Martens, Val.co
GHM GROUP – Val.co
Via Rovereto 9/11
20014 S. Ilario di Nerviano
Milano (MI) | ITALY
Phone +39 0331 53 59 20
alessandro.perego@valco.it

Netherlands
GHM Meettechniek BV
Zeeltweg 30
3755 KA Eemnes | NETHERLANDS
Phone +31 35 53805-40
Fax +31 35 53805-41
info@ghm-nl.com | www.ghm-nl.com

South Africa
GHM Messtechnik SA (Pty) Ltd
16 Olivier Street
Verwoerdpark, Alberton 1453
SOUTH AFRICA
Phone +27 74 4590040
j.grobler@ghm-sa.co.za

